

Essence

Newly designed, Huvitz continues to lead in product development combining innovation with value and performance

HNT-1/1P

Huvitz Tonometer with
Smart Puffing Control Tech



Huvitz Re:define, Re⁺create

Smoothness and Perfection HNT-1/1P

Soft & Smart Puffing, Corneal Thickness Compensation, Combined with Great Economical Value – A New Standard in Intraocular Pressure Measurement

The new tonometer HNT-1/1P measures customized intraocular pressure with smart function auto-adjustable puffing intensity.

Intuitive interface based on corneal thickness to compensate for IOP value, produces accurate measuring data immediately and effortlessly.



IOP Measurement Taking into Account Corneal Thickness – Producing Comprehensive Data



Accurate Corneal Thickness Compensation

To measure accurate intraocular pressure, simply input patient's corneal thickness on the HNT-1 to print-out compensated IOP value. (HNT-1)

When using the built-in pachymeter, available on the HNT-1P, it immediately shows compensated IOP value. (HNT-1P)

CCT(Central Cornea Thickness) Measurement

Enables accurate measuring corneal thickness by utilizing the scheidpflug method. (HNT-1P)

Visualization for Corneal Thickness Measurement

Bilateral corneal thickness by visualizing cross-section image of measured corneal thickness. (HNT-1P)

ACA(Anterior Chamber Angle) Capture

ACA cross-section capturing function helps to support the diagnosis of angle-closure which is one of the main causes of glaucoma. (HNT-1P)

Angle Measurement Function by Touch Screen

Utilize the ACA cross-section touch-screen, and the angle shows in graphic and numerical display with easy measurement. (HNT-1P)



Intuitive Gesture and Easy to Use

High Resolution 7" Color Touch-Screen

By adopting a wide color TFT LCD, it produces a vivid, high resolution (with no afterimage) image with real-time processing chip.

User-friendly and easy to use touch-screen.

More Accurate Safety Stopper Function

When pushing the safety button, it prevents contacting air nozzle to patient's eye by means of adopting auto-sensor which initializes the position of the air nozzle.

Motorized Chin Rest

User-friendly and easy to use motorized chin rest.

High Speed Internal Printer

Built-in printer, conveniently and quickly prints measured data.

Network Data Transfer Function


Send measured data to external computer by RS-232C interface cable. (EMR compatible)

Power Saving Sleep Mode Function


Automatic sleep mode when not in use.

CCT Measurement and IOP Compensation

Visualization for Corneal Thickness

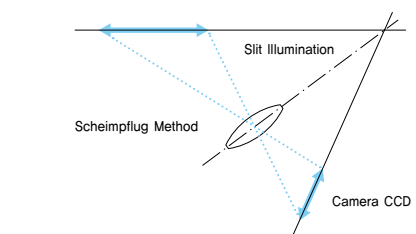

CCT(Central Cornea Thickness) Measurement Concept

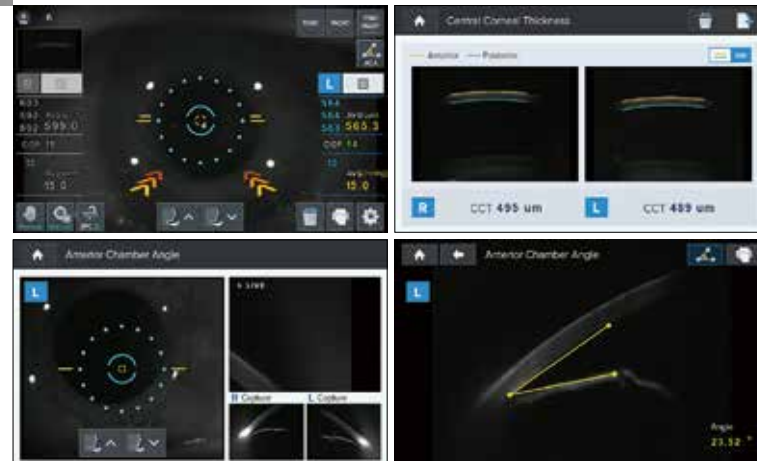

ACA Shooting

ACA Measurement



Motorized Chin Rest

Internal Printer



HNT-1/1P Huvitz Tonometer with Smart Puffing Control Tech

Specification

		HNT-1	HNT-1P
Intraocular Pressure Measurement Mode	Measurement Mode	AT3D(X, Y, Z), AT2D(X, Y), MT(Manual)	
	Measurable Range	0~60mmHg, SPC30/60mmHg	
	Measurement Value	1mmHg step (Average:0.1mmHg step)	
Cornea Thickness Measurement Mode	Pachymetry Measurement		AT3D(X,Y,Z), AT2D(X,Y), MT(Manual)
	Measurable Range		150~1300μm
	Measurement Value		1μm step
Data Memory		Measured value of ten(10) times amount for each left/right eye.	
Hardware	Built-in Printer	Thermal Line Printer	
	Power Saving Function	As stopping to measure for about 1/3/5minutes, the main power is shut, it returns as pushing button.	
	Monitor	TFT LCD Color IPS Touch Panel of 7" (800x480)	
	Electrical Power / Current	AC100~240V, 50/60Hz, 1A	
	External I/O	RS-232C (in/out)	
	Dimensions / Weight	514(W) x 262(D) x 435(H)mm / 18.5(17.5)kg	

Designs and details can be changed without prior notice for the purposes of improvement.



HUVITZ Co., Ltd. 38, Burim-ro 170beon-gil, Dongan-gu, Anyang-si, Gyeonggi-do, 14055, Republic of Korea
 Tel:+82-31-442-8868 Fax:+82-31-477-8617
<http://www.huvitz.com>

Distributed by